

ApoAlert® Caspase Colorimetric Assay Kits User Manual

Table of Contents

I. Introduction	3
II. List of Components	5
III. Additional Materials Required	5
IV. ApoAlert® Caspase-3 Assay Protocol	6
A. General Considerations	6
B. Caspase-3 Assay Procedure	7
V. ApoAlert® Caspase-8 Assay Protocol	8
A. General Considerations	8
B. Caspase-8 Assay Procedure	9
VI. Colorimetric Quantification of Protease Activity	10
VII. References	11
VIII. Related Products	12

Notice to Purchaser

This product is intended to be used for research purposes only. It is not to be used for drug or diagnostic purposes, nor is it intended for human use. Clontech products may not be resold, modified for resale, or used to manufacture commercial products without written approval of Clontech Laboratories, Inc.

Clontech, Clontech logo and all other trademarks are the property of Clontech Laboratories, Inc. Clontech is a Takara Bio Company. ©2005

I. Introduction

The **ApoAlert® Caspase Colorimetric Assay Kits** provide a simple means for detecting the activity of key caspases in mammalian cells: caspase-3 and caspase-8. These caspases are part of a large family of cysteine proteases that mediate apoptosis, or programmed cell death. Upon activation, they cleave a variety of target proteins in the cell, causing morphological and functional changes to cells undergoing apoptosis (Lazebnik et al., 1994; Casciola-Rosen et al., 1996).

Background

Apoptosis plays a fundamental role in many normal biological processes as well as in several disease states (Thompson, 1995). Apoptosis can be induced by various stimuli that all produce the same end result: systematic and orderly cell death.

Caspases cleave a variety of cellular substrates after aspartic acid residues—a characteristic that is central to their role in mammalian apoptosis. Caspases are synthesized in the cytosol of mammalian cells as inactive zymogens, which become active through intracellular caspase cascades (Cohen, 1997). The ApoAlert Caspase Colorimetric Assay Kits allow you to detect the activity of a specific caspase—caspase-3 or caspase-8—which becomes active at different stages and/or different pathways of the apoptotic process.

Caspase-3 is an active cell-death protease involved in the execution phase of apoptosis, where cells undergo morphological changes such as DNA fragmentation, chromatin condensation, and apoptotic body formation (Porter & Janicke, 1999; Zou et al., 1999). Caspase-3 is activated in response to serum withdrawal, activation of Fas, and treatment with radiation and pharmacological agents (Zou et al., 1999), as are other upstream caspases such as caspase-8 and caspase-9.

Caspase-8, the most upstream caspase in the CD95 apoptotic pathway, is activated by the signaling pathways for CD95/Fas and tissue necrosis factor (TNF; reviewed in Nicholson & Thornberry, 1997; Thornberry & Littlewood, 1998; Fraser & Evan, 1996; Cohen, 1997). Upon activation by cell-surface receptors, such as Fas, caspase-8 directly or indirectly initiates the proteolytic activities of downstream effector caspases, such as caspase-3 (Srinivasula et al., 1996; Muzio et al., 1997; Zou et al., 1999).

For a review of caspases and apoptosis, see Green & Reed, 1998.

I. Introduction continued

The ApoAlert Caspase Colorimetric Assay Kits provide a simple means for assaying caspase-3 and caspase-8 activity in mammalian cell lines. Figure 1 shows the cleavage reactions. The colorimetric assay uses the spectrophotometric detection of the chromophore p-nitroaniline (pNA) after its cleavage by caspases from the labeled caspase-specific substrates.

Inhibitors for both caspases are available for investigating the overall roles of proteases in the apoptotic process. Caspase-3 is specifically inhibited by the synthetic tetrapeptides DEVD-CHO and DEVD-fmk. DEVD-CHO is a reversible inhibitor; DEVD-fmk is irreversible. Caspase-8 is inhibited irreversibly by IETD-fmk. These inhibitors can be used as negative controls. See Related Products (Section VIII) for ordering information.

Comparing the reading of an apoptotic sample with an uninduced control allows determination of the fold-increase in protease activity. Units of protease activity can also be quantified accurately and reproducibly using a standard curve established with free pNA (Section VI).

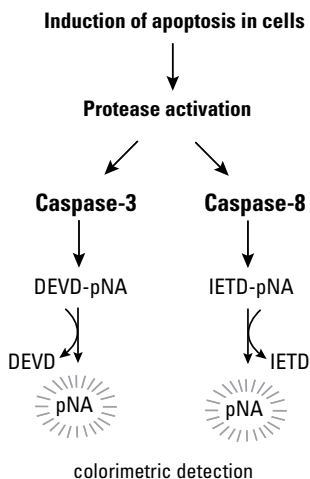


Figure 1. Detecting protease activity using the ApoAlert® Caspase Colorimetric Assay Kits. Colorimetric detection for caspase-3 and caspase-8 is performed using a spectrophotometer at 405 nm.

II. List of Components

Store Cell Lysis Buffer and 2X Reaction Buffer at 4°C after opening.

Free pNA may be stored at 4°C after first use.

Store all other components at –20°C.

Caspase-3 Colorimetric Assay Kits

25 assays (630216)	100 assays (630217)	
• 32 ml	125 ml	Cell Lysis Buffer
• 1.25 ml	5 ml	2X Reaction Buffer
• 250 µl	1 ml	DTT (1 M)
• 125 µl	500 µl	Caspase-3 Substrate, DEVD-pNA (1 mM)
• 15 µl	50 µl	Caspase-3 Inhibitor, DEVD-fmk (1 mM)
• 50 µl	200 µl	pNA (100 mM)

Caspase-8 Colorimetric Assay Kits

25 assays (630220)	200 assays (630221)	
• 25 ml	100 ml	Cell Lysis Buffer
• 2 ml	16 ml	2X Reaction Buffer
• 100 µl	800 µl	DTT (1 M)
• 25 ml	2 x 100 ml	Dilution Buffer
• 125 µl	1 ml	Caspase-8 Substrate, IETD-pNA (4 mM)

III. Additional Materials Required

The following materials are required but not supplied:

- **pNA** (Sigma #N2128) for Caspase-8 assay
- **Caspase-8 Inhibitor** (Cat. No. 630209) for Caspase-8 assay
- **Microcentrifuge** for collecting cells
- **Reaction tubes (0.5–1.5 ml) or 96-well plates**
- **Spectrophotometer and cuvettes, or microplate reader capable of measurement at 400–405 nm.**
- **DMSO**

IV. ApoAlert Caspase-3 Assay Protocol

PLEASE READ THE ENTIRE PROTOCOL BEFORE BEGINNING

A. General Considerations

- A relatively high concentration of DTT (10 mM) is required for full enzyme activity. Ensure that DTT is added to the Reaction Buffer when the assay is performed. Otherwise, unexpected low caspase activity will occur. Turbidity or particulate materials in samples can decrease assay precision.
- If not analyzed immediately, samples should be stored at -20°C or lower.
- Protect DEVD-pNA from light.
- Cover or cap all reagents when not in use, and do not mix reagents from various kits.

Reagents

- Aliquot a sufficient volume of 2X Reaction Buffer (50 μl /assay) for the number of assays you will perform. **Immediately before use**, add DTT to the aliquotted 2X Reaction Buffer to obtain a final concentration of 10 mM (1:100 dilution of stock DTT solution). **Do not** add DTT to stock 2X Reaction Buffer.
- Ensure that all buffers are completely in solution before use. After thawing, store the Cell Lysis Buffer and 2X Reaction Buffer at 4°C .

Controls

- We recommend performing three control reactions:
 - 1) a negative control on uninduced cells (Step IV.B.1)
 - 2) a control on induced cells treated with Caspase-3 Inhibitor (Step IV.B.6)
 - 3) a positive control for Caspase-3 using a known apoptosis inducer.

Cells

You may perform the assay in 0.5–1.5-ml reaction tubes or 96-well plates. Use at least 2×10^6 cells per sample. Using fewer cells may reduce the observed increase in protease activity.

IV. ApoAlert Caspase-3 Assay Protocol continued

B. Caspase-3 Assay Procedure

1. Induce apoptosis in cells using the desired method. Remember to incubate a concurrent control culture without induction.

Set up duplicate cell plates for the following samples: induced, uninduced (negative control), induced plus inhibitor (Step 6; optional), and induced without substrate (Step 8).

2. Count cells and centrifuge 2×10^6 cells at $400 \times g$ for 5 min.

Note: After removing the supernatant, you may freeze the cell pellets at -70°C and assay at a later time.

3. Resuspend cells in 50 μl of chilled Cell Lysis Buffer.

4. Incubate cells on ice for 10 min.

5. Centrifuge cell lysates in a microcentrifuge at maximum speed for 10 min at 4°C to precipitate cellular debris. Transfer the supernatants to new microcentrifuge tubes and place them on ice.

Note: At this point, samples may be assayed immediately or frozen for assay at a later time.

6. **[Optional]** In the absence of a positive control reaction with purified Caspase-3, incubate an induced sample with Caspase-3 Inhibitor before adding substrate. This reaction verifies that the signal detected by the kit is due to protease activity.

Add 50 μl of **2X Reaction Buffer/DTT Mix** (Section IV.A) and 1 μl of Caspase-3 Inhibitor (DEVD-fmk) to 50 μl of supernatant from a sample obtained in Step 5. Incubate on ice for 30 min together with the other samples. Proceed to Step 7, where you will add reaction buffer to the remaining samples.

7. Add 50 μl of **2X Reaction Buffer/DTT Mix** (Section IV.A) to each reaction. If you are using the Caspase-3 Inhibitor as a negative control, add 1 μl DMSO per 50 μl of 2X Reaction Buffer to samples without inhibitor to ensure that all samples are tested under similar conditions.
8. Add 5 μl of 1 mM Caspase-3 Substrate (DEVD-pNA; 50 μM final conc.) to each tube. Incubate at 37°C for 1 hr (or up to 3 hours maximum) in a water bath. Set up a parallel control reaction that does not contain conjugated substrate.

9. Read samples at 405 nm in a microplate reader or spectrophotometer using a 100- μl quartz cuvette. Alternatively, you may dilute the samples to 1 ml with PBS and use a standard cuvette, although doing so will weaken the signal. For quantification of protease activity, proceed to Section VI.

Note: Subtract the background reading from cell lysates and buffers from the reading of induced and uninduced samples before calculating the fold-increase in Caspase-3 activity.

V. ApoAlert Caspase-8 Assay Protocol

PLEASE READ THE ENTIRE PROTOCOL BEFORE BEGINNING

A. General Considerations

- A relatively high concentration of DTT (10 mM) is required for full activity of the enzymes. Ensure that DTT is added to the Reaction Buffer when the assay is performed. Otherwise, unexpected low caspase activity will occur. Turbidity or particulate materials in samples can decrease assay precision.
- If not analyzed immediately, samples should be stored at -20°C or lower.
- Protect IETD-pNA from light.
- Cover or cap all reagents when not in use, and do not mix reagents from various kits.

Reagents

- Aliquot a sufficient volume of 2X Reaction Buffer (50 μl /assay) for the number of assays you will perform. **Immediately before use**, add DTT to the aliquotted 2X Reaction Buffer to obtain a final concentration of 10 mM (1:100 dilution of stock DTT solution). **Do not** add DTT to stock 2X Reaction Buffer.
- Ensure that all buffers are completely in solution before use. After thawing, store the Cell Lysis Buffer and 2X Reaction Buffer at 4°C .

Controls

- We recommend performing three control reactions:
 - 1) a negative control on uninduced cells (Step V.B.1)
 - 2) a control on induced cells treated with Caspase-8 Inhibitor (630209; see Section V.B.6).
 - 3) a positive control for Caspase-8 using a known inducer of apoptosis.

Cells

You may perform the assay in 0.5–1.5-ml reaction tubes or 96-well plates. Use at least 2×10^6 cells per sample. Using fewer cells may reduce the observed increase in protease activity.

V. ApoAlert Caspase-8 Assay Protocol continued

B. Caspase-8 Assay Procedure

1. Induce apoptosis in cells using the desired method. Remember to incubate a concurrent control culture without induction.

Set up duplicate cell plates for the following samples: induced, uninduced (negative control), induced plus inhibitor (Step 6; optional), and induced without substrate (Step 8).

2. Count cells and centrifuge 2×10^6 cells at $400 \times g$ for 5 min.

Note: After removing the supernatant, you may freeze the cell pellets at -70°C and assay at a later time.

3. Resuspend cells in 50 μl of chilled Cell Lysis Buffer.

4. Incubate cells on ice for 10 min.

5. Centrifuge cell lysates in a microcentrifuge at maximum speed for 10 min at 4°C to precipitate cellular debris. Transfer the supernatants to new microcentrifuge tubes and place them on ice.

Note: At this point, samples may be assayed immediately or frozen for assay at a later time.

6. **[Optional]** In the absence of a positive control reaction with purified Caspase-8, incubate an induced sample with Caspase-8 Inhibitor (630209) before adding substrate. This reaction verifies that the signal detected by the kit is due to protease activity.

Add 50 μl of **2X Reaction Buffer/DTT Mix** (Section V.A) and 1 μl of Caspase-8 Inhibitor (IETD-fmk) to 50 μl of supernatant from a sample obtained in Step 5. Incubate on ice for 30 min together with the other samples. Proceed to Step 7, where you will add reaction buffer to the remaining samples.

7. Add 50 μl of **2X Reaction Buffer/DTT Mix** (Section V.A) to each reaction. If you are using the Caspase-8 Inhibitor as a negative control, add 1 μl DMSO per 50 μl of 2X Reaction Buffer to samples without inhibitor to ensure that all samples are tested under similar conditions.

8. Add 5 μl of 4 mM Caspase-8 Substrate (IETD-pNA; 200 μM final conc.) to each tube. Incubate at 37°C for 1 hr (or up to 3 hours maximum) in a water bath. Set up a parallel control reaction that does not contain conjugated substrate.

9. Read samples at 405 nm in a microplate reader or spectrophotometer using a 100- μl quartz cuvette. Alternatively, you may dilute the samples to 1 ml with PBS and use a standard cuvette, although doing so will weaken the signal. For quantification of protease activity, proceed to Section VI.

Note: Subtract the background reading from cell lysates and buffers from the reading of induced and uninduced samples before calculating the fold-increase in Caspase-8 activity.

VI. Colorimetric Quantification of Protease Activity

pNA Calibration Curve

1. Generate a pNA calibration curve:
 - a. Dilute 100 mM pNA solution in DMSO or in Cell Lysis Buffer to make 0, 0.5, 1, 2, and 4 mM stock solutions.
 - b. To 5 µl of each stock solution, add 95 µl of Cell Lysis Buffer to give these final concentrations:

5 µl 0 mM pNA + 95 µl of Buffer = 0 nmole pNA

5 µl 0.5 mM pNA + 95 µl of Buffer = 2.5 nmole pNA (25 µM)

5 µl 1.0 mM pNA + 95 µl of Buffer = 5.0 nmole pNA (50 µM)

5 µl 2.0 mM pNA + 95 µl of Buffer = 10 nmole pNA (100 µM)

5 µl 4.0 mM pNA + 95 µl of Buffer = 20 nmole pNA (200 µM)
 - c. Measure the five dilutions with a spectrophotometer at 405 nm and prepare a calibration curve with **x** = nmole pNA and **y** = O.D. Units (ODU).

Sample results (read at 405 nm; your results may vary):

<u>nmole pNA</u>	<u>ODU</u>
0	0.001
2.5	0.060
5.0	0.138
10.0	0.237
20.0	0.489

2. Use the slope ($\Delta\text{ODU}/\Delta\text{nmol pNA}$) of this curve to calculate units of caspase activity with the following formula:

$$\text{Units caspase activity} = \Delta\text{ODU} \times \frac{1}{\text{curve slope}}$$

ΔODU = the difference in ODU between an uninduced control and its corresponding induced sample.

VII. References

- Casciola-Rosen, L. A., Nicholson, D. W., Chong, K. R., Rowan, K. R., Thornberry, N. A., Miller, D. K. & Rosen, A. (1996) Apopain/CPP32 cleaves proteins that are essential for cellular repair: a fundamental principle of apoptotic cell death. *J. Exp. Med.* **183**:1957–1964.
- Cohen, G. M. (1997) Caspases: the executioners of apoptosis. *Biochem. J.* **326**:1–16.
- Fraser, A., & Evan, G. (1996) A license to kill. *Cell* **85**:781–784.
- Green, D. R. & Reed, J. C. (1998) Mitochondria and apoptosis. *Science* **281**:1309–1312.
- Lazebnik, Y. A., Kaufmann, S. H., Desnoyers, S., Poirier, G. G. & Earnshaw, W. C. (1994) Cleavage of poly(ADP-ribose) polymerase by a proteinase with properties like ICE. *Nature* **371**:346–347.
- Muzio, M., Salvesen, G. S., & Dixit, V. M. (1997) FLICE induced apoptosis in a cell-free system. *J. Biol. Chem.* **272**:2952–2956.
- Nicholson, D. W. (1996) ICE/CED3-like proteases as therapeutic targets for the control of inappropriate apoptosis. *Nature Biotechnol.* **14**:297–301.
- Nicholson, D. W., & Thornberry, N. A. (1997) Caspases: killer proteases. *Trends Biochem. Sci.* **22**:299–306.
- Porter, A. G. & Janicke, R. U. (1999) Emerging roles of caspase-3 in apoptosis. *Cell Death Diff.* **6**:99–104.
- Srinivasula, S. M., Ahmad, M., Fernandes-Alnemri, T., Litwack, G., & Alnemri, E. S. (1996) Molecular ordering of the Fas-apoptotic pathway: the Fas/APO-1 protease MCH-5 is a CrmA-inhibitable protease that activates multiple Ced-3/ICE-like cysteine proteases. *Proc. Natl. Acad. Sci. USA* **93**:14486–14491.
- Thompson, C. B. (1995) Apoptosis in the pathogenesis and treatment of disease. *Science* **267**:1456–1462.
- Thornberry, N.A. & Littlewood, Y. (1998) Caspases: Enemies Within. *Science* **281**:1312–1316.
- Zou, H., Li, Y., Liu, X. & Wang, X. (1999) An APAF-1 cytochrome c multimeric complex is a functional apoptosome that activates procaspase-9. *J. Biol. Chem.* **274**:11549–11556.

VIII. Related Products

For a complete listing of all Clontech products,
please visit www.clontech.com

Apoptosis Inducing Reagent

- Human TNF- α 630203

ApoAlert® Apoptosis Inhibiting Reagents

- Caspase-3 Inhibitor, DEVD-CHO 630204
- Caspase-3 Inhibitor, DEVD-fmk 630207
- Caspase-8 Inhibitor, IETD-fmk 630209
- Caspase Inhibitor, VAD-fmk 630208

ApoAlert® Apoptosis Assay Kits

- Caspase 3 Fluorescent Assay Kit 630214
630215
- Caspase 8 Fluorescent Assay Kit 630218
630219
- Mitochondrial Membrane Sensor Kit 630106
- DNA Fragmentation Assay Kit 630107
630108
- Annexin V-FITC Apoptosis Kit 630109
630110
- Cell Fractionation Kit 630105
- Caspase-9/6 Fluorescent Assay Kit 630211
630212
- Caspase-3 Assay Plate 630223
630224
- Caspase Profiling Plate 630225
630226

ApoAlert® Apoptosis Vectors

- pDsRed2-Bid 632419

Notes

Notes

Notes
